# Vibraluxe PRO®

 *Instruction Manual*

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 *Manual de instrucciones*

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

 *Mode d'emploi*

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International patents pending. All imitations will be prosecuted. Vibraluxe PRO © is a registered EU/CTM trade mark.

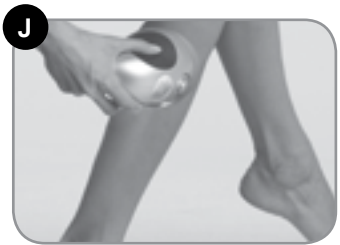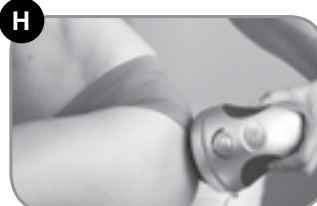Thank you for purchasing Vibraluxe Pro®. Vibraluxe Pro® is based on professional appliances designed to help reduce the appearance of cellulite, flabbiness and water retention. Vibraluxe Pro® provides an intense, deep massage of the tissues below the surface of the skin using both circular rotation and vibration to help improve circulation, make the skin look smoother, tone muscles, shape contours, help reduce the appearance of cellulite and eliminate excess water. Its ergonomic design makes it your perfect partner for frequent use.

**IMPORTANT: PLEASE READ THESE INSTRUCTIONS CAREFULLY BEFORE USING THIS APPLIANCE**

#### **WARNING:**

- Do not use Vibraluxe Pro® if you are pregnant or if you have a heart condition or physical injury.
- **IMPORTANT:** Do not use the appliance if you have implants such as a pacemaker, an IUD, surgical screws, etc..
- Keep Vibraluxe Pro® out of reach of children.
- Do not use Vibraluxe Pro® if you have thrombosis or are being treated for any medical condition.
- Do not use Vibraluxe Pro® for 1 hour before or after meals.
- If you have had a serious illness, consult your doctor before using Vibraluxe Pro®.
- Do not use Vibraluxe Pro® near water or with wet hands.
- Stop use immediately if you notice any form of discomfort.
- Unplug after use.
- Do not dismantle, repair or modify

#### **Vibraluxe Pro®.**

- Vibraluxe Pro® is not a medical or therapeutic device.
- In case of doubt, consult your doctor about using the appliance.
- Do not use the product on areas where there is inflammation or irritation of the skin.
- Do not use the appliance for more than 5 minutes on the same area of the body.
- Do not use the product prior to going to sleep, since the massage may have a stimulating effect.
- Do not use the appliance while lying down on the bed.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- Do not use the product without one of the heads.
- The device cannot be used if the head is broken.
- When using the product without the mesh

protective cover, we recommend tying your hair back to prevent it from getting tangled in the rotating head and causing physical injury.

Using Vibraluxe Pro® has many beneficial effects on your body, thanks to its combination of rotation and vibration which produces a deep massage of the tissues below the surface of the skin:

- It shapes contours and tones muscles: its rotating motion and vibration reactivates your muscles. You will achieve a toned and slimmer look with Vibraluxe Pro®, which includes also a diet plan for enhanced results.

- It reduces the appearance of cellulite: it produces a deep, intense massage. The skin looks smoother with a reduction of the orange-peel effect.

- It helps fight water retention: it improves circulation and encourages the elimination of liquids and toxins.

- Helps relax areas with muscle tension: using your Vibraluxe Pro® on the low intensity levels provides a feeling of wellbeing in areas with muscle tension.

The set includes (A):

- Vibraluxe Pro® Massager
- Smooth Head: for a relaxing massage
- Large-Rippled Head: to massage, tone and shape large areas
- Small-Rippled Head: to massage, tone and shape small areas
- Multi-Node Head: for a toning, firming massage
- Pumice Head: to smooth rough skin on feet
- Protective Cover: to fit over the head attachment for a gentler massage

- Travel case
- Diet plan
- Instruction manual.

#### RECOMMENDATIONS FOR USE:

- If desired, firming or anti-cellulite creams can be used with your Vibraluxe Pro®.
- Do not use your Vibraluxe Pro® for more than five minutes on the same area.

#### Accessories (B):

- Smooth Head (1): perfect for a relaxing massage to soothe accumulated tension.
- Large-Rippled Head (2): perfect for a toning, firming massage of stubborn areas where fat often accumulates such as tummy, thighs and buttocks.
- Small-Rippled Head (3): perfect for a toning, firming massage and for treating localized fat on small areas like arms and hips.
- Multi-Node Head (4): perfect for a toning massage and a firming effect.
- Pumice Head (5): perfect for smoothing rough skin on feet. Do not use creams with this attachment.
- Cover (6): fit the cover over the head attachment for massaging very sensitive areas and making massage easier when Vibraluxe Pro is used directly on the skin. Only use the cover with the smooth or rippled heads. It is recommended that you use the cover when massaging the shoulder area to avoid hair becoming tangled in the massager. Do not use creams when the cover is on.

#### ADJUSTABLE INTENSITY LEVEL

The adjustable intensity dial (C) lets you adapt the massage to the area of the body you are treating. Low intensity is recommended for more sensitive areas like arms, calves, tummy, legs with varicose veins, back and shoulders. High intensity is recommended for areas with localized fat like thighs and buttocks.

#### FREQUENCY OF USE

Initially you should use your Vibraluxe Pro® once a day, then 2-3 times a week to maintain results. Use for a maximum of 5 minutes on each area.

#### INSTRUCTIONS FOR USE:

- Plug the appliance into an electrical socket.

- Fit the head that you wish to use by pushing it onto the appliance. If desired, attach the protective cover by pushing it onto the head. To remove the protective screen, first align the tab on the cover with the triangle, then turn the cover and remove. (D)

- Select the massage intensity setting on the intensity dial.
- When the massage is over, turn the power dial to OFF.
- Unplug the appliance.

#### RECOMMENDED USES

- E- Tummy: perfect for shaping the tummy area and helping reduce localized fat.
- F- Shoulders and Back: perfect for soothing accumulated tension and relaxing the back area.
- G- Upper Arms: for toning the underside of your arms.
- H- Buttocks and Hips: perfect for toning and helping to reduce stubborn fat.
- I- Thighs: perfect for shaping, toning and helping reduce localized fat.
- J- Calves: perfect for soothing tight muscles.
- K- Feet: perfect for smoothing rough skin on feet.

#### CLEANING, MAINTENANCE AND STORAGE:

After unplugging the appliance, clean with a soft, damp cloth. Do not use soap. Keep the appliance in a cool, dry place, preferably stored in the travel case with the head attachments.

#### QUALITY WARRANTY

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country.

This warranty does not cover damages resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.

#### Australia only

This product has the benefit of certain Consumer Guarantees prescribed by the Australian Consumer Law which provide protection for consumers. There is no express warranty for this product in Australia. The above paragraph refers to other countries.

#### DISPOSAL

With the aim of protecting the environment, if you decide to dispose of this device, the battery or the current adapter, use the special containers for electrical and electronic devices.

#### CERTIFICATES AND REGULATIONS

This product complies with Low Voltage Directive 2014/35/EU, Electromagnetic Compatibility Directive 2014/30/EU, RoHS Directive 2011/65/EU and its amendment Directive (EU) 2015/863 and LFGB and REACH Certifications.

#### TECHNICAL SPECIFICATIONS

Power: 25 W  
Frequency: 50 Hz  
Voltage: 220V - 240V  
Minimum rpm: 1000  
Maximum rpm: 3000

This appliance has a double insulating system, so it does not need ground connection. Please make sure that the tension of the mains corresponds to the value indicated on the appliance's information badge.

#### ATTENTION DISPOSAL OF MATERIALS

A symbol of a crossed-out wheeled bin means you should find out about and follow local regulations about disposing of this kind of product.

 Do not dispose of this product as you would other household waste. Dispose of this device in accordance with the corresponding local regulations.

Electrical and electronic devices contain hazardous substances that can have harmful effects on the environment and/or human health and should be recycled properly.

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ISL GmbH
Best Direct (International) Benelux B.V.
Best Direct International Ltd,
ISL ITALY S.r.l.