GYM
form®

TOTAL FITNESS

Instruction Manual

Manual de instrucciones

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Mode d'emploi

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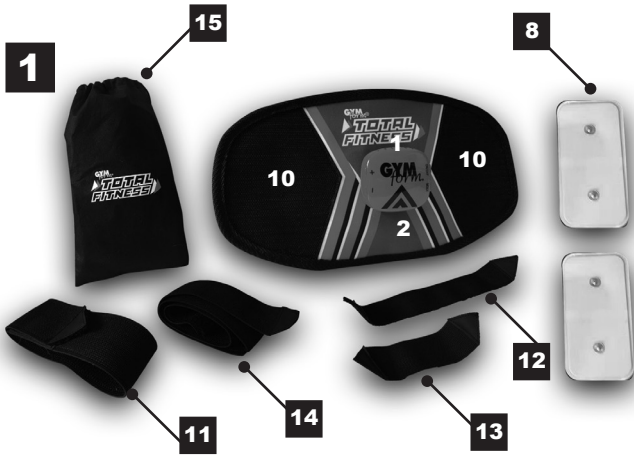
Congratulations on your purchase of the Gymform® TOTAL FITNESS. Your safety is our primary concern, so please ensure that you familiarize yourself with all the information in this manual before using your Gymform® TOTAL FITNESS. Anyone who intends to use this device must read and understand the instructions provided. An informed user is a safe user. Thank you for purchasing the Gymform® TOTAL FITNESS. We hope you enjoy your sessions.

WARNINGS

Before embarking on any exercise programme, including the Gymform® TOTAL FITNESS, consult your general practitioner and your physiotherapist. Ensure that you read the instructions, warnings and precautions carefully before using your Gymform® TOTAL FITNESS.

The Gymform® TOTAL FITNESS is designed for healthy adults. However, the device should be used with care. Read the following guidelines and instruction manual before using the product.

- Applying electrodes in the thorax area may increase the risk of cardiac fibrillation.

- Using the Gymform® TOTAL FITNESS while connected to high-frequency surgical equipment may cause burns at the contact points of the Gymform® TOTAL FITNESS and may result in damage to the product.

- Operation of a short-wave therapeutic device or microwave oven in the general vicinity (e.g. 1 metre away) may cause instability in the output power of the Gymform® TOTAL FITNESS.

- This is a powerful electrical muscle stimulation device and may cause injury if you start with the most advanced exercise programme right away. Start with a lower level and gradually

increase the intensity.

- Do not use on children or animals.
- Do not use while driving or sleeping.
- Do not use in the bathroom or in other places with a high amount of moisture or humidity.
- The batteries should be removed from the device before it is disposed of.
- The batteries should be safely disposed of.

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

- Batteries that are not rechargeable should not be recharged.

- Dead batteries should be removed from the device and safely disposed of.

- Remove the batteries from the device if you do not plan to use it for an extended period of time.

- Do not short-circuit the battery poles.

DO NOT use the Gymform® TOTAL FITNESS:

- If you have an electronic implant (such as a pacemaker) or suffer from another heart problem.
- If you are, or think you might be, pregnant.
- If you have cancer, epilepsy or are

under medical supervision for a cognitive disorder.

- If you are breastfeeding.
- If you have a high fever.

Before using Gymform® TOTAL FITNESS, women should wait:

- At least six months after the birth of your child (after consulting your doctor).
- A month after having an intrauterine device (IUD) fitted.
- At least three months after having a caesarean section (after consulting your doctor).
- Until after the first days of your period.

Before using Gymform® TOTAL FITNESS, you should obtain permission from your doctor or physiotherapist:

- If you suffer from back problems. In this case, the device should be used only at low intensity levels.
- If you suffer from a serious disease or injury not mentioned in this section.
- If you were operated on recently.
- If you are diabetic and take insulin.
- Do not place the electrodes on your neck, chest or back in such a way that electricity passes through your heart (e.g. one electrode on your chest and one on your back), genitals, scars, moles, varicose veins or head.
- Do not place on recent scars, inflamed skin, infected areas or areas that are susceptible to acne, thrombosis or other vascular problems (e.g. varicose veins), or any part of your body with limited sensitivity or hypersensitivity.
- Do not place on areas with injuries or limited mobility (e.g. fractures or sprains).
- Do not place the product directly on

metal implants. Place on the closest muscle.

Possible Adverse Reactions:

There have been isolated reports of skin reactions caused by the use of electro-stimulation devices, including allergies, prolonged skin redness and acne. If any of these symptoms occur, switch off the device and consult your doctor.

If you have never exercised before or have not exercised for a long time: Consult your general practitioner.

Precautions

- Do not use Gymform® TOTAL FITNESS if you are drowsy or while sleeping.
- Do not use Gymform® TOTAL FITNESS if your body is in contact with water.
- Keep the device out of reach of children.
- Do not use Gymform® TOTAL FITNESS in combination with any other device that transmits electrical current to your body or with any other machine.
- Do not connect the metal snap fasteners or electrodes to any other object.
- Do not modify the device.
- Use only the electrodes supplied by Industex.
- Do not attach two electrodes together.
- If any of the electrodes are not attached properly, you will notice an itching sensation.

The Set Contains

- 1 Gymform® TOTAL FITNESS console with six programmes
- 1 Console belt

- 1 Adjustable elastic strap to wrap around the waist, abdominals, chest and buttocks (75 cm)
- 1 Elastic strap extension for large waists (38 cm)
- 1 Extra elastic strap for thighs (30 cm)
- 1 Extra elastic strap for calves, legs and arms (20 cm)
- 1 Instruction manual
- 1 Travel bag
- 1 Screwdriver

PRODUCT DESCRIPTION (Picture 1):

1. Gymform® TOTAL FITNESS unit
2. Console belt
3. ON/OFF button
4. + Button (to increase intensity level)
5. - Button (to reduce intensity level)
6. Mode button (programme selection)
7. Programme indicator lights
8. Intensity level indicator
9. Electrodes
10. Adjustable fasteners
11. Strap for waist, abdominals, chest and buttocks (75 cm)
12. Strap extension for large waists (38 cm)
13. Strap for thighs (30 cm)
14. Strap for calves, legs and arms (20 cm)
15. Travel bag

Assembly

Open the battery compartment cover (batteries not included) on the back of the unit (1) (see Picture 2). Insert two CR2032 batteries. Use the screwdriver to remove the battery compartment cover: insert it into the screw head and turn until the screw comes out. After inserting the batteries, replace the cover and position the screw. Then use the screwdriver to tighten it. Insert both batteries in accordance with the polarity marks. The polarity is indicated on the battery compartment for proper placement.

Putting on the Belt

Attach the Gymform® TOTAL FITNESS (1) to the console belt (2) with the two snap fasteners. Ensure that they are firmly attached. If not, the belt will not work properly.

Attach the two electrodes (9) to the inner part of the belt. Use the snap fasteners to firmly attach the electrodes.

Remove the protective plastic backing from the belt electrodes. Do not throw the protective plastic backing away because you will need it later.

See Picture 3 to decide on the area of your body you want to work on. Choose the right elastic strap for the body part you want to work on.

- Adjustable elastic strap (11) to wrap around the waist, abdominals, chest and buttocks (75 cm)
- Strap extension for large waists (12) (38 cm)
- Extra elastic strap (13) for thighs (30 cm)
- Extra elastic strap (14) for calves, legs and arms (20 cm)

Attach each end of the elastic strap to the console belt using the adjustable fasteners.

Note: The strap extension (12) is only necessary for using the belt on large waists. Attach one end of the strap extension (12) to the elastic strap (11) and the other end to the console belt (2) using the adjustable fasteners.

When you have finished using the belt, remove the electrodes one at a time to ensure they do not stick to the skin. Pull on them from the edge and then replace the protective plastic backing on the electrodes. This will help extend the useful life of the electrodes (9).

Programme Selection

Press the ON/OFF button (3) to turn on the unit (1). When Gymform® TOTAL FITNESS is first switched on, it starts at Programme 1 with an intensity of 0. Use the MODE button (6) to select other programmes. Press the MODE button (6) once to change to Programme 2, once again for Programme 3, and so on. After Programme 6, Programme 1 will be displayed again. It's easy to know which programme is selected by looking at the programme indicator lights (7).

Programme	Function
1	Massage
2	Aerobic Programme
3	Muscle Toning
4	General Workout
5	Anti-Cellulite
6	Drainage

Adjusting the Intensity Level

To reduce the intensity level, briefly and repeatedly press the - Button (5). To increase the intensity level, briefly and repeatedly press the + Button (4). There are 10 intensity levels to choose from. It's easy to know which intensity level is selected by looking at the intensity level indicator lights (8). If this is your first electro-stimulation session, we

recommend that you start at Level 1 for the first 10 minutes and then gradually increase the intensity as you become more familiar with the device. When you have finished the session, press the ON/OFF button (3) to switch off the unit (1). We recommend using the product for 10 minutes a day in a single area.

Automatic Power Off

For your safety, the automatic period for each programme is 10 minutes. Gymform® TOTAL FITNESS will automatically switch itself off 10 minutes after being switched on, regardless of the level or programme selected.

If you want to stop the session before that time, press the ON/OFF button (3) to turn off the unit (1).

Automatic Reset

When you finish the programme, the device will switch itself off automatically. When you switch it on again, you will have to programme it again.

WE RECOMMEND THAT YOU REPLACE THE ELECTRODES AFTER EVERY 30 USES. IF YOU NOTICE THE ELECTRODES HAVE PREMATURELY LOST THEIR ADHESIVE POWER, WET THEM SLIGHTLY WITH A DAMP SPONGE.

NOTE: If you feel heat in the areas where you have attached the electrodes or if the contractions become weak during the session, it is likely that the batteries need to be replaced. Remove them and insert new ones. If this does not solve the problem, the electrodes have probably reached the end of their useful life. Switch off the Gymform® TOTAL FITNESS, take it off, fit new electrodes to the belt and then resume your session. For your convenience, Gymform® TOTAL FITNESS will automatically set the intensity level to 0 if you change programmes.

TECHNICAL SPECIFICATIONS

Temperature: 0~+35°C

<85% RH

700-1060 hPa

Continuously operating device.

Power supply: 2 x 3 V CR2032 batteries (not included)

Voltage: 2 x 3 V (6 V total)

Programme	Impulse Width (µs)	Frequency (Hz)	Voltage (V)
1	100	2	36 V
2	100	90	36 V
3	100	120	36 V
4	100	10	36 V
6	100	50	36 V

Maximum current of the electrodes supplied: 1.397 mA (at 4000 Ohm).

Dimensions: 92 x 82 x 15 mm (main unit).

Weight: 76 g (main unit without batteries).

This symbol indicates that the device is protected against electrostatic discharge by BF-type double insulation.

ISO 7000-1641 User Guide

Values measured at 4 kΩ

Safety classification based on the type of protection against electric shock; internally powered device with a type BF applied part.

Picture 1
User's Guide

CLEANING AND MAINTENANCE

DO NOT WASH THE GYMFORM® TOTAL FITNESS UNIT: WIPE WITH A DRY OR DAMP CLOTH.

Remove the batteries from the Gymform® TOTAL FITNESS if you do not plan to use it for an extended period of time.

CLEANING AND MAINTENANCE OF ADHESIVE ELECTRODES

-Do not immerse in water.

-Keep in a cool, dry place at room temperature.

-Do not remove the insulating label.

-After using the device, always replace the protective plastic backings on the adhesive electrodes.

- Store in the bag supplied.
- Clean by wiping with wet fingers.
- Do not remove the insulating label.

PRECAUTIONS WITH THE BATTERIES

Ensure the batteries are properly installed before switching on the unit.

Turn off the unit to change the batteries.

Remove the batteries if you are not going to use the device for an extended period of time. Take the necessary precautions when disposing of the batteries. Do not expose them to fire. Made in China

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Para su seguridad, el periodo automático de cada programa es de 10 minutos. Gymform © TOTAL FITNESS se apagará automáticamente 10 minutos después de haberlo encendido, independientemente del Nivel y Programa configurados.
ISL GmbH
ISL ITALY S.r.l.
Best Direct (International) Benelux B.V.
Best Direct International Ltd,